



stretcher series 5110 general transport stretcher

The full function 5110 Series general transport stretcher provides a versatile yet economical patient transport solution. This is a non-hydraulic adjustable stretcher.



RAIL SYSTEM

Easy to wipe clean and disinfect, full length tuck-away side rails allow for easier patient transfer by creating a zero transfer gap. Rails feature a release latch centrally located on either side of the litter top under the patient surface.



FOWLER ASSIST

Pneumatic assist fowler provides additional assistance for the safe and easy raising and lowering of the patient.

Pneumatic Assist Fowler

400 N (90#) force



ZERO TRANSFER GAP

Tuck-Away side rails allows the caregiver to position the 5110 Series Stretchers flush with the bed for a zero transfer gap and easier patient transfer. No expensive or cumbersome transfer devices are required.



HYDRAULIC PEDALS

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Color-coded pedals and function graphics are quickly and easily understood.



COLUMN SUPPORT

Sturdy steel column design with five fixed litter top height positions to choose from: 28", 29", 30", 31" or 32".

PATIENT LIFTS

The 5400 Series Stretcher has 57 $\frac{7}{8}$ "H" of clearance between floor and cowlings to easily accommodate and utilize patient lifts.



TRUE DIRECTION STEERING

For improved stretcher control engage the True Direction Steering system. The steer caster located at the patient's right foot end of the stretcher will lock once the steering pedal is applied.




STANDARD FEATURES

- 500 lbs. weight capacity
- True Direction Steering (TDS)
- Side rails - full length, single release
- Center column patient support
- Non-hydraulic
- 5 fixed litter top height options: 28" to 32"
- Single pneumatic cylinder assists in placing the litter top in Trendelenburg/ Reverse Trendelenburg position.
- Pneumatic assist fowler 400 N (90#) force
- Brake & Steer activated foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- High impact corner bumpers
- 8" easy mobility casters
- Latex-free

SPECIFICATIONS

Model Number	5110 ^N	5110 ^W
Weight Capacity	500 lbs. (227 kg)	500 lbs. (227 kg)
Overall Length	83" (211 cm)	83" (211 cm)
Overall Width	29.5" (75 cm)	33.5" (85 cm)
Fixed Height Range		
Low	28" (71 cm)	28" (71 cm)
	29" (73 cm)	29" (73cm)
	30" (76 cm)	30" (76 cm)
	31" (79 cm)	31" (79 cm)
	32" (81 cm)	32" (81 cm)
Litter Positioning		
Fowler Backrest	0°-90°	0°-90°
Trendelenburg	+12°	+12°
Reverse Trendelenburg	-7°	-7°
Patient Surface	77" x 27" (196 x 69 cm)	77" x 31" (196 x 77 cm)
Side Rails	15" x 63 ⁵ / ₈ " (38 x 162 cm)	15" x 63 ⁵ / ₈ " (38 x 162 cm)
Casters	8" (20 cm)	8" (20 cm)
Warranty	1 year parts and labor	1 year parts and labor


OPTIONS

Accessory	N	W	Accessory	N	W	Accessory	N	W
X-Ray Board			IV Pole Options:			Additional Options:		
Full.....	750203	750204	Fixed Height.....	2101	same	Bed Exit Monitor.....	750216	same
Partial.....	750205	750206	Adjustable / Removable.....	5984001	same	Storage Bag Hook.....	750215	same
Mattresses:			Adjustable / Attached.....	750501	750502	Paper Roll Holder.....	5966001	5966003
3" Standard (black).....	5828002	5831002	Fold Down.....	59175001	same	Quick Release O ₂ Holder.....	750213	same
3" Premium (blue).....	5850002	5852002	IV Pole Caddy.....	59163001	same	O ₂ Tank Holder (Vertical).....	5981001	same
4" Standard (black).....	5829002	5832002	Push/Pull Bar Options:			Head/Foot Extension.....	5996001	same
4" Premium (blue).....	5851002	5853002	Push/Pull Handles (CuVerro®).....	750223	750224	Head/Foot 3" Pad.....	5708001	same
4" Pressure Relief.....	5833002	5834002	Push/Pull Bar (Removable).....	5985001	5986001	Stirrups Set.....	59100001	same
5" Premium (blue)†.....	5855002	5856002	Bar Pad.....	5704001	5705001			
† <i>Special order only.</i>			Chart Holder.....	5967001	same			
Restraint Straps Options:			Monitor Tray Board.....	59115001				
Velcro (pair).....	5903001	same	59116001					
Buckle (pair).....	5904001	same	(includes chart holder 5967001)					
Side Rail Options:			Knee Flex					
Rail Pad (pair).....	5716002	same	Hand Crank.....	750201	750202			
			Hydraulic.....	750208	750209			
								
			Leg Lift					
				Knee Flex				
					Trauma			



IV CADDY

Quickly and easily secure an IV pole to a stretcher using Pedigo's IV Caddy. Allows one caregiver to manage the safe transport of the patient, stretcher and IV equipment. Stores under litter top when not in use.

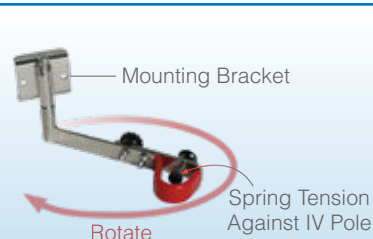

FOLD DOWN IV POLE