



stretcher series 5400 universal procedure stretcher

The 5400 Series universal transport stretcher is versatile, yet economical. It's the right stretcher for ED, OR, PACU, or any department where cost, efficiency and quality are required.



RAIL SYSTEM

Easy to wipe clean and disinfect, full length tuck-away side rails allow for easier patient transfer by creating a zero transfer gap. Rails feature a release latch centrally located on either side of the litter top under the patient surface.



FOWLER ASSIST

Pneumatic assist fowler provides additional assistance for the safe and easy raising and lowering of the patient.

Pneumatic Assist Fowler

400 N (90#) force



ZERO TRANSFER GAP

Tuck-Away side rails allows the caregiver to position the 5110 Series Stretchers flush with the bed for a zero transfer gap and easier patient transfer. No expensive or cumbersome transfer devices are required.



HYDRAULIC PEDALS

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Color-coded pedals and function graphics are quickly and easily understood.




COLUMN SUPPORT

Steel center column support with single hydraulic lift cylinder design uses sixteen precision ball-bearings to smoothly raise and lower the patient.

PATIENT LIFTS

The 5400 Series Stretcher has 5 $\frac{7}{8}$ "H" of clearance between floor and cowling to easily accommodate and utilize patient lifts.

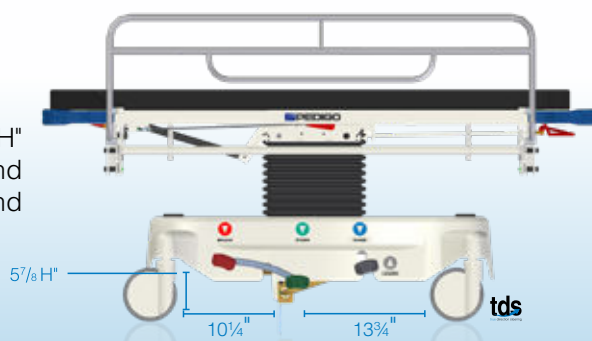

TRUE DIRECTION STEERING

For improved stretcher control engage the True Direction Steering system. The steer caster located at the patient's right foot end of the stretcher will lock once the steering pedal is applied.




STANDARD FEATURES

- 500 lbs. weight capacity
- True Direction Steering (TDS)
- Side rails - full length, single release
- Center column design with precision ball bearings
- Single Trendelenburg assist cylinder (positive/negative - pneumatic)
- Pneumatic assist fowler 400 N (90#) force
- Dual side-mounted hydraulic foot pedals
- Brake & Steer activated foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- High impact corner bumpers
- 8" easy mobility casters
- Latex-free

SPECIFICATIONS


Model Number	5400 ^N	5400 ^W
Weight Capacity	500 lbs. (227 kg)	500 lbs. (227 kg)
Overall Length	83" (211 cm)	83" (211 cm)
Overall Width	29.5" (75 cm)	33.5" (85 cm)
Height Range		
Low	23" (58 cm)	23" (58 cm)
High	32.75" (83 cm)	32.75" (83 cm)
Litter Positioning		
Fowler Backrest	0°-90°	0°-90°
Trendelenburg	+12°	+12°
Reverse Trendelenburg	-7°	-7°
Patient Surface	77" x 27" (196 x 69 cm)	77" x 31" (196 x 77 cm)
Side Rails	15" x 63 ⁵ / ₈ " (38 x 162 cm)	15" x 63 ⁵ / ₈ " (38 x 162 cm)
Casters	8" (20 cm)	8" (20 cm)
Warranty	1 year parts and labor	1 year parts and labor

OPTIONS

Accessory	N	W	Accessory	N	W	Accessory	N	W
X-Ray Board			IV Pole Options:			Additional Options:		
Full	750203	750204	Fixed Height	2101	same	Bed Exit Monitor	750216	same
Partial	750205	750206	Adjustable / Removable ..	5984001	same	Extra Release Pedal - Hydraulic for stretchers with TDS	750105	same
Mattresses:			Adjustable / Attached.....	750501	750502	iS6™ Instant Steer 6 Wheel	750103	same
3" Standard (black).....	5828002	5831002	Fold Down	59175001	same	Extra Release Pedal - Hydraulic for stretchers with iS6™ ..	750106	same
3" Premium (blue).....	5850002	5852002	IV Pole Caddy	59163001	same	Storage Bag Hook	750215	same
4" Standard (black).....	5829002	5832002	Push/Pull Bar Options:			Paper Roll Holder	5966001	5966003
4" Premium (blue).....	5851002	5853002	Push/Pull Handles (<i>CuVerro®</i>)	750223	750224	Quick Release O ₂ Holder	750213	same
4" Pressure Relief	5833002	5834002	Push/Pull Bar (Removable) ..	5985001	5986001	O ₂ Tank Holder (Vertical) ..	5981001	same
5" Premium (blue)†	5855002	5856002	Bar Pad	5704001	5705001	Head/Foot Extension	5996001	same
† <i>Special order only.</i>			Chart Holder	5967001	same	Head/Foot 3" Pad	5708001	same
Restraint Straps Options:			Monitor Tray Board	59115001	59116001	Stirrups Set	59100001	same
Velcro (pair)	5903001	same	<i>(includes chart holder 5967001)</i>					
Buckle (pair)	5904001	same	Knee Flex					
Side Rail Options:			Hand Crank	750201	750202			
Rail Pad (pair)	5716002	same	Hydraulic	750208	750209			
			  					
			Leg Lift Knee Flex Trauma					



SIDE RAIL PADS

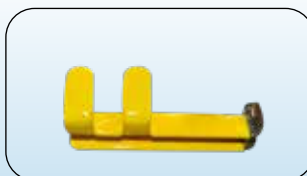

STORAGE BAG HOOK



RESTRAINT STRAPS