




stretcher series  
**7500<sup>NE</sup>**  
 eye | neck | head stretcher



The new 7500<sup>NE</sup> Stretcher is designed specifically for eye, neck and head surgery. Standard features include: center column support, retractable push handles, Quick-Release O<sub>2</sub> holder, and **IS6** (instant steer six wheel steering).

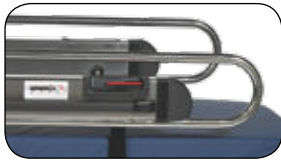

**RAIL SYSTEM**

Easy to wipe clean and disinfect in the ongoing battle with germs and bacteria. This 56" stainless steel side rail is infinitely adjustable and can be locked at any height without interfering with the surgical process or patient entry and exit from the stretcher.



**O<sub>2</sub> TANK HOLDER**

Pedigo's ergonomic Quick-Release O<sub>2</sub> Tank Holder is side mounted and conveniently located under the patient surface to reduce caregiver bending.

(US Patent 7,865,984 B2)




**SIX WHEEL STEERING**

The two Instant Steer wheels are centered directly under the chassis for increased control and maneuverability. Dual ball bearing wheels feature a softer wheel tread for improved grip and traction.



**HYDRAULIC PEDALS**

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Color-coded pedals and function graphics are quickly and easily understood.

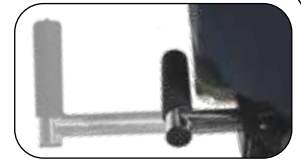

**PUSH/PULL HANDLES**


Easily extend or retract push handles by engaging the safety release lock and rotate each handle upward.



**head rest articulation**




Head Rest Button Function

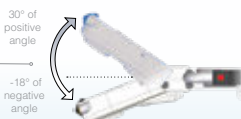

Head Rest Angular Articulation



1.5" Up

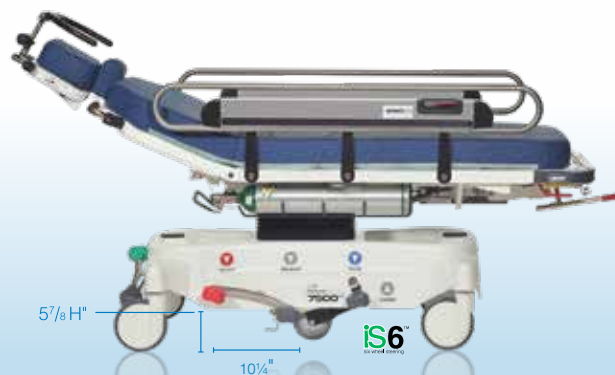

**PATIENT LIFTS**

The 7500<sup>NE</sup> Series Stretcher has 57<sup>7</sup>/<sub>8</sub>"H of clearance between floor and cowling to easily accommodate and utilize patient lifts.

**INSTANT STEER SIX WHEEL STEERING**

The two **iS6<sup>TM</sup>** wheels are centered directly under the chassis for increased maneuverability and control. Dual ball bearing wheels feature a softer wheel grip and traction.

## STANDARD FEATURES

- 750 lbs. weight capacity
- Stainless steel glide side rails
- Pneumatic articulating head and neck rests
- Pads - 4" Head and 6" Neck
- 4" Premium mattress (part # 750719)
- Pneumatic adjustable fowler
- Integrated push/pull handles with latch release
- Four IV pole receptacles
- Four restraint strap positions
- Quick-Release O<sub>2</sub> Tank Holder (US Patent 7,865,984 B2)
- Hydraulic recovery knee flex
- Hydraulic height adjustment - dual side-mounted foot pedals
- **IS6**™ Instant Steer - Six Wheel Steering
- Single column hydraulic cylinder design
- Brake foot pedals
- 8" easy mobility casters
- Latex-free

## SPECIFICATIONS

|                           |                          |
|---------------------------|--------------------------|
| <b>Model Number</b>       | 7500 <sup>NE</sup>       |
| <b>Weight Capacity</b>    | 750 lbs. (340 kg)        |
| <b>Overall Length</b>     | 88" (224 cm)             |
| <b>Overall Width</b>      | 34" (86 cm)              |
| <b>Height Range</b>       |                          |
| Low                       | 23" (59 cm)              |
| High                      | 32.75" (83 cm)           |
| <b>Litter Positioning</b> |                          |
| Fowler Backrest           | 0°-90°                   |
| Trendelenburg             | +12°                     |
| Reverse Trendelenburg     | -7°                      |
| <b>Head Rest Range</b>    |                          |
| Up                        | 1.5"                     |
| Down                      | 1.5"                     |
| Total Range               | 3"                       |
| <b>Head Rest Tilt</b>     |                          |
| Range                     | +30° -18°                |
| <b>Patient Surface</b>    | 82" x 27" (208 x 69 cm)  |
| <b>Side Rails</b>         | 12.5"x 56" (32 x 142 cm) |
| <b>Casters</b>            | 8"(20 cm)                |
| <b>Warranty</b>           | 1 year parts and labor   |

## OPTIONS

### Accessory

#### Electric Fowler Accessories:

|                                  |        |
|----------------------------------|--------|
| Battery Charging Station .....   | 750711 |
| Battery for 7500-NE Fowler ..... | 750720 |

#### Fowler Option:

|                             |                          |
|-----------------------------|--------------------------|
| Electric Back Support ..... | 750705                   |
|                             | <i>Factory Installed</i> |

#### Head Rest Extension.....750704

#### IV Pole:

|                                    |          |
|------------------------------------|----------|
| Fixed Height .....                 | 2101     |
| Height Adjustable (Tethered) ..... | 750702   |
| Fold Down .....                    | 59175001 |

#### Miscellaneous:

|                                   |        |
|-----------------------------------|--------|
| Extended Trend Release Handle ..  | 750708 |
| Patient Belongings Bag Hook ..... | 750706 |
| Raised Brake Pedals .....         | 750721 |

### Accessory

|                                    |          |
|------------------------------------|----------|
| Transfer Pad .....                 | 750724   |
| Extra Hydraulic Release Pedal..... | 750106   |
| Stainless Accessories Bar.....     | 75070006 |

#### Pad Kit - Supplemental .....750722

|                            |  |
|----------------------------|--|
| Contoured Head Pad         |  |
| Head Pad, 2"H x 8.5" x 9"  |  |
| Head Pad, 3"H x 8.5" x 9"  |  |
| Head Pad, 4"H x 5" x 3.25" |  |
| Neck Pad, 5"H x 5" x 3.25" |  |

#### Pad Kit - Complete .....750701

|                           |  |
|---------------------------|--|
| Contoured Head Pad        |  |
| Head Pad, 2"H x 8.5" x 9" |  |
| Head Pad, 3"H x 8.5" x 9" |  |
| Head Pad, 4"H x 8.5" x 9" |  |

### Accessory

|                            |
|----------------------------|
| Neck Pad, 4"H x 5" x 3.25" |
| Neck Pad, 5"H x 5" x 3.25" |
| Neck Pad, 6"H x 5" x 3.25" |

#### Push/Pull Bar:

|  |          |
|--|----------|
| Push/Pull Bar ( <i>Removable</i> ) .....   | 5985001  |
| Bar Pad .....                              | 5704001  |
| Chart Holder.....                          | 5967001  |
| Monitor Tray Board .....                   | 59115001 |
| (includes chart holder 5967001)            |          |
| O <sub>2</sub> Tank Holder (Vertical)..... | 5981001  |

#### Restraint Straps:

|                     |         |
|---------------------|---------|
| Velcro (pair) ..... | 5903001 |
| Buckle (pair) ..... | 5904001 |

#### Surgeon Wrist Rest:

|                               |          |
|-------------------------------|----------|
| Non-swiveling Wrist Rest..... | 750703   |
| Swiveling Wrist Rest .....    | 750709   |
| 3" Height Extension .....     | 750709-1 |



HEAD REST EXTENSION



ELECTRIC FOWLER OPTION



HEAD & NECK PADS